


20 August 2021

# COVID-19 TESTING

## Specimen requirement: One Nasopharyngeal swab

### NASOPHARYNGEAL SWAB COLLECTION

- Wear PPE: mask, gloves, eye protection, gown
- Rest back of patient's head against wall or bed
- Insert swab 4-5 cm into nose in horizontal direction
- Rotate swab 3 x each way

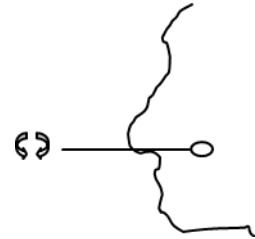

## COVID-19 swab packs

- Each pack contains one swab and one tube
- Swab is **flexible**
- Tube contains **virus transport medium** (VTM or VTF)

*Actual swabs/tubes may vary slightly from those pictured*



- There is no requirement to double bag swabs for COVID-19 testing
- COVID-19 swabs **must be bagged separately from other laboratory specimens** such as bloods and other swabs, and have a separate laboratory request form
- COVID-19 swabs will routinely be tested for SARS-CoV-2 virus only
- If additional respiratory swab examinations are required, these must be ordered at time of initial request. There will be no add on testing.
- Also accepted: Sputum (but this requires aerosol precautions for collection)